



SPECIAL NEEDS STROLLERS CATALOGUE





The MOUSE special needs stroller is designed for the smallest children who are unable to walk or sit in a correct position as a result of injuries, diseases or motor dysfunctions. The stroller is suitable for children at a height of 80-110 cm and a weight of 30 kg.

Thanks to the backrest tilt and recline feature as well as adjustable torso and head pads you may easily adjust MOUSE to the required position and ensure the user with optimal protection and

The shock-absorbing suspension on all four, large, gel-filled wheels together with the tilt-foot facilitate travel on uneven terrain. Modern and stylish design is an additional advantage of

INDICATIONS:

paralysis and paresis of lower limbs of various origins paralysis and paresis of the trunk weakening of postural muscles responsible for maintaining the correct body position cerebral palsy myelomeningocele muscular dystrophies



colour options.

















TILT IN SPACE SEAT

simple mechanism



8°-23°

4 positions

RECLINING BACKREST

wide range of adjustments



16°-88°

4 positions

ADJUSTABLE AND ANGLE TILT FOOTREST

simple and quick adjustment



84°- 173°

6 positions

SEAT MINIMIZING KIT

reducing insert



The SEAT MINIMIZING KIT in MOUSE is a comfortable, one-peace insert that allows you to reduce both width and depth of the seat.

STABILIZATION SYSTEM

adjustable pads



The head and torso support consists of two pairs of independent pads that are height, width and angle adjustable. Its wide range of adjustments and simplicity gives you a fast and effective adaptation to the user's individual needs.

LYING POSITION

wide range of backrest recline



REVERSIBLE SEAT

arrange your travel easily



The design of the stroller allows you to easily rotate the seat opposite to the driving direction. Thanks to this feature you can keep an eye on your child all the time

FOLDING LIGHTWEIGHT FRAME


ease of storage, transportation and dismantling



The construction of main elements allows you to dismantle Mouse into 3 separate units (frame, seat and wheels).

DRIVING SUPPORT

comfort for the bumpiest journeys

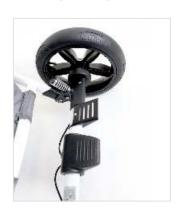

The stroller is equipped with a shock-absorbing suspension on all four wheels and an adjustable rear suspension in the rear wheels, swivel front castor wheels with option to lock, central brake and a tilt-foot.





Specifications

TECHNICAL DATA OF A STROLLER READY TO DRIVE



seat depth 22*/26 cm

ubackrest height 52 cm

footrest adjustable within 14,5 - 23 cm

- **total weight 17,2 23,2 kg**
- **b** total leng**t**h 95 124 cm
- total height 122 134 cm
- total width 60 cm

dimension after folding: frame L84 x W 58 x H 25 cm / seat L72 x W 42 x H 39 cm

USER PARAMETERS



maximum weight 30 kg





The MEWA special needs stroller is designed for children who are unable to walk or sit in a correct position as a result of injuries, diseases or motor dysfunctions. The stroller is suitable for children at a height of 90-125 cm and a weight of 40 kg.

Thanks to the backrest tilt and recline feature as well as adjustable torso and head pads you may easily adjust MEWA to the required position and ensure the user with optimal protection and stability.

The shock-absorbing suspension on all four, large, gel-filled wheels together with the tilt-foot facilitate travel on uneven terrain. Modern and stylish design is an additional advantage of this stroller.

INDICATIONS:

paralysis and paresis of lower limbs of various origins paralysis and paresis of the trunk weakening of postural muscles responsible for maintaining the correct body position cerebral palsy myelomeningocele muscular dystrophies

colour options.



FRAME

spray painted / glossy





TILT IN SPACE SEAT

simple mechanism



5,5°-21,5°

4 positions

SEAT MINIMIZING KIT

reducing inserts



The SEAT MINIMIZING KIT in MEWA is a set of 4 comfortable, independent inserts that allows you to reduce both width and depth of the seat.

RECLINING BACKREST

wide range of adjustments



8°-86°

4 positions

STABILIZATION SYSTEM

adjustable pads



The head and torso support consists of two pairs of independent pads that are height, width and angle adjustable. Its wide range of adjustments and simplicity gives you a fast and effective adaptation to the user's individual needs.

ADJUSTABLE AND ANGLE TILT FOOTREST

simple and quick adjustment



64°- 180°

5 positions

LYING POSITION

wide range of backrest recline



REVERSIBLE SEAT

arrange your travel easily



The design of the stroller allows you to easily rotate the seat opposite to the driving direction. Thanks to this feature you can keep an eye on your child all the time

FOLDING LIGHTWEIGHT FRAME

ease of storage, transportation and dismantling



The construction of main elements allows you to dismantle MEWA into 3 separate units (frame, seat and wheels).


DRIVING SUPPORT

comfort for the bumpiest journeys



The stroller is equipped with a shock-absorbing suspension on all four wheels and an adjustable rear suspension in the rear wheels, swivel front castor wheels with option to lock, central brake and a tilt-foot.



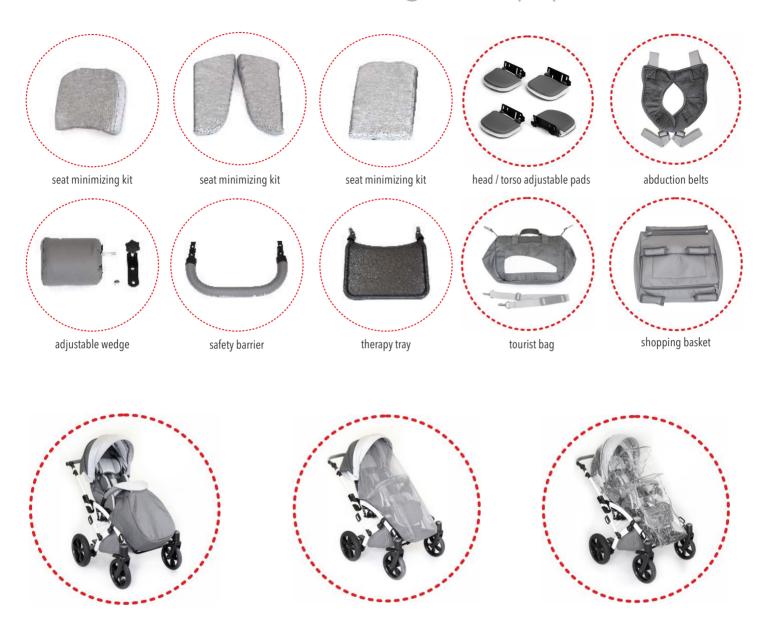

Specifications

TECHNICAL DATA OF A STROLLER READY TO DRIVE



winter sleeping bag

seat depth 24*/28 cm

ubackrest height 62 cm

footrest adjustable within 17 - 32 cm

total weight 20 - 28 kg

mosquito net

- **b** total leng**t**h 94 126 cm
- total height 124 129 cm
- total width 65 cm

dimension after folding: frame L 84 x W 64 x H 26 cm / seat L 88 x W 43 x H 42 cm

USER PARAMETERS

🙀 growth 90 - 125 cm

rain cover

maximum weight 40 kg





The GRIZZLY special needs stroller is designed for the children who are unable to walk or sit in a correct position as a result of injuries, diseases or motor dysfunctions. The stroller is suitable for children at a height of 85-140 cm and a weight of 45 kg.

Grizzly is equipped with the "I GROW WITH YOU" system that allows for a smooth extension of the backrest height and seat depth as your child grows. Moreover Grizzly offers you a comprehensive stabilization system consisting of adjustable pads at the level of the head, hips and chest – together with a fully adjustable footrest, tilt in space seat and reclining backrest we bring you a stroller that is tailored to fit your child's needs in 100% and will serve for years.

Grizzly is crash tested and it ensures safety while transporting children in specially adapted vehicles.

INDICATIONS:

paralysis and paresis of lower limbs of various origins paralysis and paresis of the trunk weakening of postural muscles responsible for maintaining the correct body position cerebral palsy myelomeningocele muscular dystrophies

colour options.





UPHOLSTERY


synthetic





spray painted / glossy



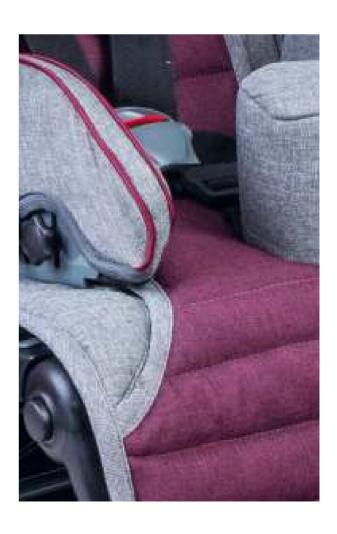

TILT IN SPACE SEAT simple mechanism



0°-30° 4 positions

I GROW WITH YOU

extensive adjustment options



"I GROW WITH YOU" system allows to extend: - the seat depth up to 37 cm

- backrest height up to 70 cm (no tools required) Designed to adapt and evolve!

RECLINING BACKREST

wide range of adjustments



0°-81° 6 positions

STABILIZATION SYSTEM

safety and stability



A complete support system consisting of a four independent pairs of pads (head, torso, hips and lateral) guarantees maximum of adjustment possibilities as well as highest comfort and safety during travel.

ADJUSTABLE AND ANGLE TILT FOOTREST

simple and quick adjustment



55°- 180° 5 positions

LYING POSITION

wide range of backrest recline



REVERSIBLE SEAT

arrange your travel easily



The design of the stroller allows you to easily rotate the seat opposite to the driving direction. Thanks to this feature you can keep an eye on your child all the time

FOLDING LIGHTWEIGHT FRAME

ease of storage, transportation and dismantling



The construction of main elements allows you to dismantle GRIZZLY into 3 separate units (frame, seat and wheels).

DRIVING SUPPORT

comfort for the bumpiest journeys



The stroller is equipped with a shock-absorbing suspension on all four wheels and an adjustable rear suspension in the rear wheels, swivel front castor wheels with option to lock, central brake and a tilt-foot





Specifications

TECHNICAL DATA OF A STROLLER READY TO DRIVE



seat depth 25 - 37 cm




footrest adjustable within 23 - 40,5 cm







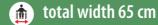

dimension after folding: frame L 80 x W 64 x H 35 cm / seat L 72 x W 43 x H 47 cm

USER PARAMETERS



growth 85 - 140 cm



maximum weight 45 kg

* with adjustable pads





The Yeti special needs stroller is designed for larger/older children who are unable to walk or sit in a correct position as a result of injuries, diseases or motor dysfunctions. The stroller is suitable for users at a height of 130-175 cm and a weight of 80 kg.

Its heavy-duty construction is equipped with the "I GROW WITH YOU" system (smooth extension of the backrest height and the seat depth as your child grows) and adjustable tension system which allows to adjust the proper upholstery tension both on the seat and the backrest guarantees you a support for many years.

Moreover, Yeti's solid suspension along with a fast folding mechanism, tilted seat, reclining backrest and a wide range of accessories significantly simplify the everyday use irrespective of the weather conditions and time of the year. Modern and practical design is his additional asset.

INDICATIONS:

paralysis and paresis of lower limbs of various origins paralysis and paresis of the trunk weakening of postural muscles responsible for maintaining the correct body position cerebral palsy myelomeningocele muscular dystrophies

colour options.





FRAME spray painted / matt













TILT IN SPACE SEAT

simple mechanism

RECLINING BACKREST

wide range of adjustments

ADJUSTABLE AND ANGLE TILT FOOTREST

simple and quick adjustment



9,5°-29°

5 positions



3°- 99°

8 positions



80°- 180°

6 positions

SEAT ADJUSTMENTS

reducer cushions



Reducer cushions set consisting of 4 independent velcro pads offers an additional adaptation of the seat interior through simple attaching of the pads at desired position. Easy adjustment of seat width and depth at your fingertips.

I GROW WITH YOU

extensive adjustment options



- "I GROW WITH YOU" system allows to extend:
- the seat depth up to 45 cm
- backrest height up to 85 cm (tools required) Designed to adapt and evolve!

LYING POSITION

wide range of backrest recline



WINDSTOPPER

multifunctional upholstery



YETI is equipped with a multifunction windstopper that not only comes in handy during windy conditions but also protects against unnecessary arm movement and limits visual stimulation.

FOLDING LIGHTWEIGHT FRAME

ease of storage and transportation



The one-piece construction of YETI stroller allows for a quick and easy folding which facilitates daily use, storage and transport.

DRIVING SUPPORT

comfort for the bumpiest journeys



The stroller is equipped with adjustable rear suspension in the rear wheels, swivel front castor wheels with option to lock, central brake and a tilt-foot





Specifications

TECHNICAL DATA OF A STROLLER READY TO DRIVE



seat depth 31 - 45 cm

backrest height 65 - 85 cm

footrest adjustable within 28,5 - 43,5 cm

total weight 24 - 28,6 kg



total leng**t**h 114,5 - 146 cm



total height 117 - 138 cm



total width 66 cm



USER PARAMETERS



growth 130 - 175 cm



maximum weight 80 kg

* with reducer cushions



REHAB TECHNOLOGY DISTRIBUTION

Rehab Technology Distribution (RTD)

W: www.rtdaus.com.au

E: therapy@rtdaus.com.au

P: (02) 9838 8869